


Product Catalogue Australia & New Zealand 2024


Enteral Nutrition Feeding Tubes, Accessories & Gravity Systems



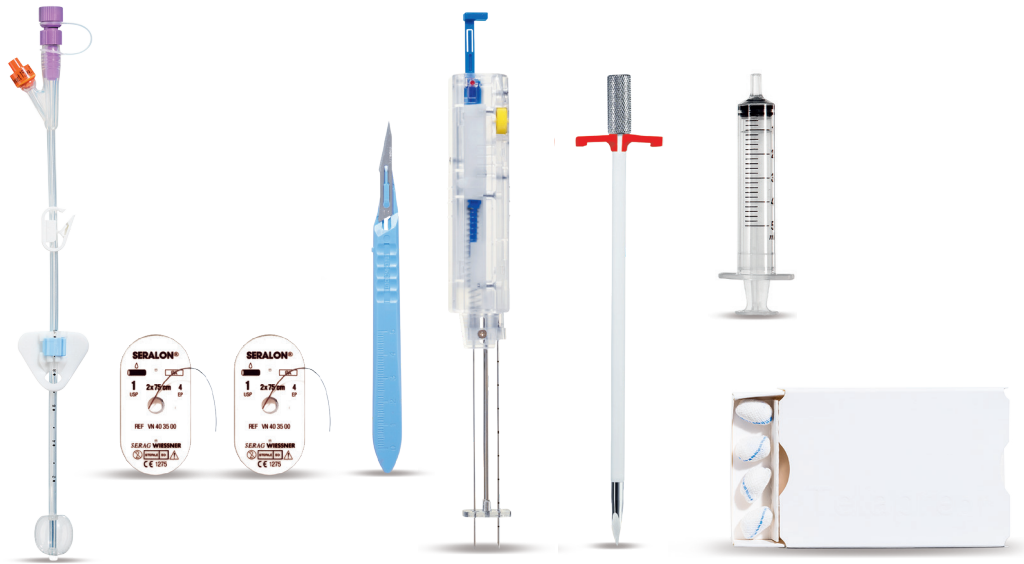

Freka® Pexact II, FR 15, ENFit

Initial placement tube

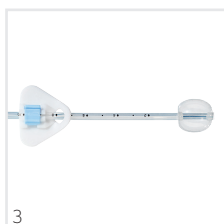

gastric

Direct Puncture Gastrostomy with secure and reliable gastropexy

- Enteral connector ENFit for a secure and tight connection
- For long-term intra-gastric feeding if direct puncture of the stomach under endoscopic control (stenosis endoscope) is required
- Gastropexy Device and suture material for secure fixing of the stomach to the internal abdominal wall (Fig. 1)
- Scalpel for point incision
- Three-corner Trocar with peel-away sheath FR 16 (Fig. 2)
- Balloon tube FR 15 of silicone rubber with ENFit Connector (Fig. 3)
- 5 ml Luer syringe for filling balloon (5 ml sterile H₂O)
- Easy change to Freka GastroTube four weeks after placement possible within a fully formed and healed stoma

Gastropexy with innovative suturing method

- Improved fixation device enables single handed operation
- Operation of the device while monitoring the endoscopic image
- Secure fixation of stomach to abdominal wall*
- Extra corporal suture technique
- Reduced risk of peritoneal infection by direct access*
- Short, straight stoma tracts

*Naoki Okumura et al., Gastroenterol Rep (Oxf). 2015 Feb; 3(1):69-71

Percutaneous endoscopic gastrostomy with Funada-style gastropexy greatly reduces the risk of peristomal infection

*Toh Yoon E W et al. PEG using a novel large-caliber introducer, Endoscopy International Open 2016

Luer-free
Latex-free
DEHP-free*

*not intentionally made with DEHP

	Connector	External ø	Length	Sales unit	Art.-no.
Freka Pexact II, FR 15, ENFit	ENF	5.0 mm	23 cm	1	7601365
ENF = ENFit					